

USER GUIDE

BEHIND-THE-EAR HEARING INSTRUMENTS



ReSound



rediscover hearing

Hearing instrument type designations for models included in this user guide are: M70-80, FCC ID: X26M70-80 and M60, FCC ID: X26M60. Please see page 10 for list of models referring to both types.

Statement:

This device complies with Part 15 of the FCC Rules.

Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from the one in which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Changes or modifications can void the user's authority to operate the equipment

Intended use

Generic air-conduction hearing instruments are wearable sound-amplifying devices intended to compensate for impaired hearing. The fundamental operating principle of hearing instruments is to receive, amplify, and transfer sound to the ear drum of a hearing impaired person.

List of countries:

Products without wireless functionality are intended for worldwide sales.

Products with wireless functionality are intended for sale in countries within the European Economic Area as well as Switzerland.

Specification of restrictions: You are not allowed to operate the equipment within 20 km of the centre of Ny Ålesund, Norway.

The products are in compliance with the following regulatory requirements:

- In EU: the device conforms to the Essential Requirements according to Annex I of Council Directive 93/42/EEC for medical devices (MDD) and essential requirements and other relevant provisions of Directive 1999/5/EC (R&TTE). The declaration of conformity may be consulted at www.resound.com
- In US: FCC CFR 47 Part 15, subpart C, section 15.249.
- Other identified applicable international regulatory requirements in countries outside the EU and US. Please refer to local country requirements for these areas.
- Products are categorized as receiver category 2 according to EN 300 440

Introduction

Congratulations on the purchase of your new hearing instruments. ReSound's innovative sound technology and design, combined with the customized programming selected by your hearing care professional, will make hearing a more enjoyable experience. Hearing instruments will enable you to hear sounds that you may not have heard in years because of your hearing loss. Practice and a positive attitude are important in learning to use hearing instruments. Your ReSound instruments have been adjusted according to your individual hearing loss and needs. Some people adjust quickly to wearing hearing instruments in their ears and hearing new sounds; other people may need more time.

Please read this manual carefully in order to wholly benefit from the use of your hearing instruments. With proper care, maintenance, and usage, your hearing instruments will aid you in better communication for many years. Ask your hearing care professional if you have any questions.

Hearing instrument model: _____

☐ Model 67: Battery size 312

☐ Model 77: Battery size 13

☐ Model 88: Battery size 13

Tube type: _____

Dome size: _____

Left serial number: _____

Right serial number: _____

Contents

Statement	2
Intended use	3
List of countries	3
Introduction	4
Description of the hearing instruments	8
Getting started	10
On/Off function	10
SmartStart	10
Inserting/Replacing the battery	10
Low battery indicator	11
Sports lock	12
Inserting/removing hearing instruments	12
Operation of the hearing instrument	14
Volume control (optional)	14
Programme button (optional)	15
Flight mode	17
Telephone use	19
Listen to radio or TV	19
Cellular phones	19
PhoneNow	20
Telecoil (optional)	22
Tele-loop systems	22
Direct audio input (optional)	22

Connecting/Disconnecting audio boots	23
Care and maintenance	24
Daily maintenance	25
Cleaning earmoulds	25
Cleaning thin tubes and domes	25
How to apply domes	26
General precautions	27
General warnings	27
For ReSound Alera TS with	
Tinnitus Sound Generator	29
Important Notice to Prospective Users	29
Prescription Use for Sound Generator Users	29
CAUTION – Tinnitus Sound Generator	30
General warnings/or precautions	30
Technical specifications	30
Available sounds	31
Battery warning Information	31
Troubleshooting guide	32
Technical Data	36
Warranty and repairs	38
Temperature test, transport and storage information	38

Hearing instruments of type M70-80 are available in the following variants:

AL 987-DVIW

AL 977-DVIW Open - AL 977-DVIW

AL 787-DVIW

AL 777-DVIW Open - AL 777-DVIW

AL 587-DVIW

AL 577-DVIW Open - AL 577-DVIW

AL 487-DVIW

AL 477-DVIW Open - AL 477-DVIW

ALT 987-DVIW

ALT 977-DVIW Open - ALT 977-DVIW

ALT 787-DVIW

ALT 777-DVIW Open - ALT 777-DVIW

ALT 587-DVIW

ALT 577-DVIW Open - ALT 577-DVIW

Hearing instruments of type M70-80e are available in the following variants:

VO988-DW

VO977-DW Open, VO977-DW

VO788-DW

VO777-DW Open, VO777-DW

Hearing instruments of type M60 are available in the following variants:

AL 967-DIW Open - AL 967-DIW

AL 767-DIW Open - AL 767-DIW

AL 567-DIW Open - AL 567-DIW

AL 467-DIW Open - AL 467-DIW


ALT 967-DIW Open - ALT 967-DIW

ALT 767-DIW Open - ALT 767-DIW

VO967-DW Open, VO967-DW

VO767-DW Open, VO767-DW


1. Programme button
2. Battery compartment & On/Off switch
3. Sound outlet
4. Front microphone inlet
5. Back microphone inlet
6. Volume control (optional)
7. Direct audio input
8. Left/Right indicator (Left=Blue/Right=Red)
9. Battery lock (optional)
10. Model
11. Manufacturer
12. Serial number
13. Sports lock

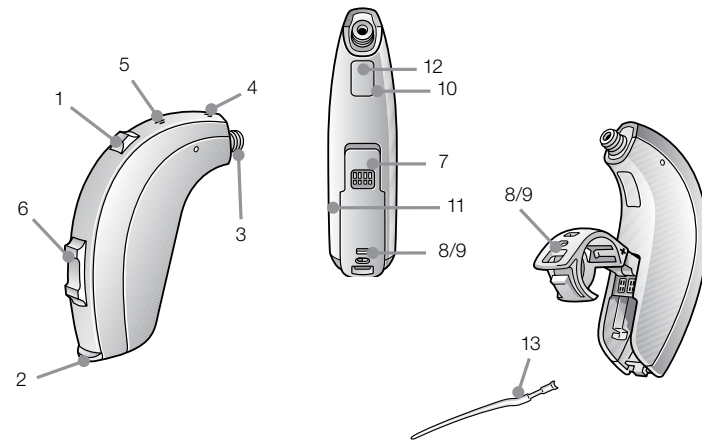

Getting started

On/Off function

1. When the battery door is closed, the hearing instrument turns on, and the default program will be activated.
2. To turn off the hearing instrument, open the battery door. Many individuals can use their fingernail to pull it open.

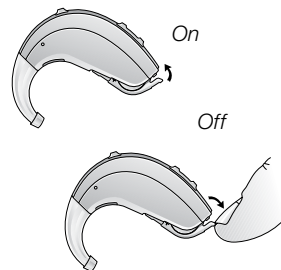


SmartStart

Hearing instruments can be turned on once you have placed them on your ears. If you prefer to turn them on just prior to placing them on your ear, your hearing care professional can activate a function called SmartStart. This function will delay the time in which the hearing instruments turn on by several seconds after the battery compartment is closed. With SmartStart, a beep will be heard for each second of the delay period.

Inserting/Replacing the battery

1. Open the battery door completely by using your fingernail.
2. Remove the used battery if present. Insert the new battery with the positive side in the correct position. Please refer to page 7 for information on appropriate battery type/size for your hearing instrument.
3. Gently close the battery door.

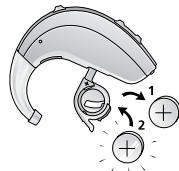

Tip:

1. Always use new Zinc-Air batteries that have a minimum remaining shelf life of 1 year.
2. Whenever the hearing instruments are not in use, remember to turn them off to avoid unnecessary battery consumption.

Low battery indicator

Your hearing care professional can activate a low battery indicator in your hearing instruments. The hearing instrument will reduce amplification and play a melody if battery power gets too low. This signal will recur every five minutes until the hearing instrument automatically switches off. The timing of the low battery indicator can vary slightly, depending on the type of battery used. It is recommended that you keep spare batteries on hand.

Low battery indicator (instruments paired with accessories only)

Active usage of the ReSound Unite accessories (Remote Control, Phone Clip, TV streamer and Mini Microphone) requires more battery power from the hearing instruments than when they are working without accessories which means that battery life is highly dependent on the amount of wireless accessory usage. When the battery in the hearing instrument has depleted to a level at which use of ReSound Unite TV Streamer, Phone Clip and Mini Microphone cannot be supported, the hearing instrument will play two sets of descending tones. After this, your hearing instrument and ReSound Unite Remote Control will continue to work as usual, but you will not be able to use your ReSound TV Streamer, Phone Clip, and Mini Microphone. At some point the battery level will also be too low to support the remote control as well and you will once again hear the descending tones. The hearing instruments will continue to work as usual. Once a new battery is inserted, full operation of the accessories will resume.

Sports lock


The Sports lock will be applied or adjusted by your hearing care professional.


Inserting/Removing hearing instruments


Insertion (custom earmoulds)

1. Hold the earmould between your thumb and index finger and position its sound outlet in your ear canal.
2. Slide the earmould all the way into your ear with a gentle, twisting movement.
3. Turn the top part of the earmould gently backwards and forwards so that it tucks behind the fold of skin above your ear canal.
4. Move the earmould up and down and gently press to ensure it is positioned correctly in the ear. Opening and closing your mouth can ease insertion.
5. Make sure the hearing instrument is seated behind the ear.

By experimenting, an easier method may be discovered. With proper insertion, hearing instruments should fit snugly but comfortably. If the hearing instruments cause irritation of the ears, contact your hearing care professional.


 **Note:** Never attempt to modify the shape of the hearing instrument, earmoulds, or tubing yourself.

 **Tip:** It may be helpful to pull your ear up and outward with your opposite hand during insertion.

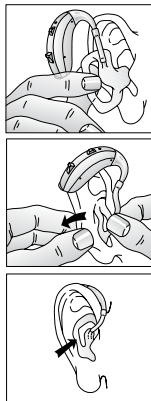


Insertion (domes/custom earmoulds with thin tubes)

1. Hold the thin tube where it bends, and gently push the dome into the ear canal. Push the dome far enough into the ear canal so that the thin tube lies flush with the head.
2. It is important that the tube and the dome fit correctly into your ear.
3. When the dome is placed correctly, you should not be able to see the thin tube sticking out when facing a mirror.

 **Note:** You should never attempt to bend or modify the shape of the thin tube.


Removal (custom earmoulds)

1. Hold a portion of your earmould towards the back of the ear.
2. Pull the earmould outward and simultaneously rotate the earmould forward.
3. Consult your hearing care professional if you have difficulties removing the hearing instrument.

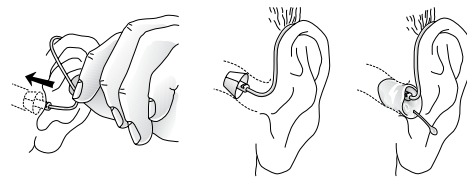

Removal (domes/custom earmoulds with thin tubes)

- 1. Hold the thin tube with your thumb and forefinger and pull the tube outward.
- 2. For thin tube custom earmoulds, grasp the removal string and pull the earmould outward.
- 3. Consult your hearing care professional if you have difficulties removing the hearing instrument.

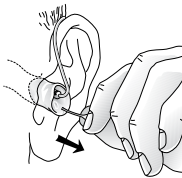


Operation of the hearing instrument

Volume control (optional)

The volume control will allow the volume of hearing instruments to be increased or decreased.

- 1. To increase the volume, push the volume control up.
- 2. To decrease the volume, push the volume control down.


When volume is increased or decreased, a beep signal will be heard for each incremental change. When the upper or lower limits of the volume range are reached, a beep signal with a longer duration will be heard.

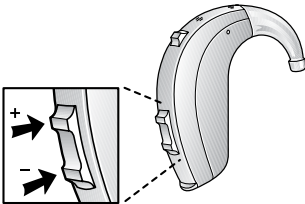

If you have two hearing instruments with the synchronization function enabled, volume control adjustments to one instrument will automatically repeat in the second instrument. When a volume control adjustment is made in one instrument, you will hear a confirmation beep. A beep in the second instrument will follow.

Programme button

Depending on your experience level with hearing instruments, individual hearing needs, and the type of listening environments you experience, your hearing care professional may activate additional programmes. If additional programmes have been activated, the following list explains how they work.

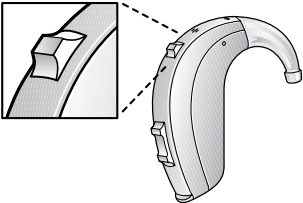

Programme	Description of when to use
1	
2	
3	
4	

- 1. You can switch between programmes by pushing the push button once.
- 2. You will then hear one or more beeps. The number of beeps indicates which programme you have selected (one beep = programme one, two beeps = programme two, etc.).
- 3. When the hearing instruments are turned off and then back on, the hearing instrument always returns to the default setting (programme one).

If you have two hearing instruments with the synchronization function enabled, program changes to one instrument will automatically repeat in the second instrument. When a program change is made in one instrument, you will hear the same amount of confirmation beeps in the second instrument.



Flight mode*

When boarding a flight **or entering an area where RF transmitters are prohibited**, wireless functionality must be deactivated, as it is not allowed to radiate radio signals during flights or in otherwise restricted areas.

For Alera wireless hearing instruments follow the following steps to enter and leave flight mode:

It is possible to disable wireless operation by opening and closing the battery compartment of the hearing instrument while at the same time pressing the push button.

When disabled manually, wireless operation may be re-enabled by opening and closing the battery compartment normally, (i.e. without at the same time pressing the push button).

For Verso wireless hearing instruments follow the following steps to enter and leave flight mode:

It is possible to disable wireless operation by performing the following sequence of operations. Note that this operation is the same for both hearing instruments with and without a program button.

1. Close the battery door (HI is turned on)
2. Open the battery door within 10 seconds of doing operation #1 (HI is turned off)


**Wireless models only*

3. Close the battery door (HI is turned on 2nd time)
4. Open the battery door within 10 seconds after doing operation #3 (WI is turned off 2nd time)
5. Close the battery door (HI is turned on 3rd time)

Telephone use

Finding the optimal position for holding a telephone may require practice for some individuals, and one or more of the following suggestions may be helpful.

1. Hold the telephone as you would normally.
2. Hold the telephone towards the top of the ear (closer to where the microphones are).
3. If whistling occurs, it may take a few seconds of holding the telephone in the same position before the hearing instrument eliminates the feedback.
4. Any whistling may also be decreased by holding the telephone slightly away from the ear.
5. Depending on your individual needs, your hearing care professional may activate a programme specifically for telephone use.

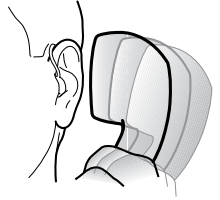

Listen to radio or TV

When listening to the TV or the radio, start out by listening to news commentators since they usually speak clearly, then try other programmes. If you find it difficult to listen to TV or radio, your hearing care professional will be able to give you advice on available accessories to enhance your listening capabilities for TV and radio.

Cellular phones

Your hearing instrument is designed to comply with the most stringent Standards of International Electro-magnetic Compatibility. However, not all cell phones are hearing instrument compatible. The varying degree of disturbance can be due to the nature of your particular cellular phone or of your wireless telephone service provider.

If you find it difficult to obtain a good result while using your cellular phone, your hearing care professional will be able to give you advice on available accessories to enhance listening capabilities.


PhoneNow

The PhoneNow function, allows your hearing instrument to automatically switch to your telephone programme when a telephone receiver is raised to the ear. When the telephone receiver is removed from the ear, the hearing instrument automatically returns to the previous listening programme.

Placement of PhoneNow magnets

Place PhoneNow magnet on your telephone receiver to allow operation of the PhoneNow function. In order to place PhoneNow magnet properly:

1. Clean the telephone receiver thoroughly.
2. Hold the telephone vertically, in a position similar to when making a telephone call.
3. Place the magnets just below the telephone receiver. Make sure not to cover the microphone openings. If necessary, move the magnet to another position to improve ease of use and comfort while speaking.
4. If you are not satisfied with the strength of PhoneNow, you can reposition the PhoneNow magnet or add additional PhoneNow magnets.

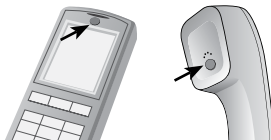

PhoneNow usage

Telephones can be used in a normal manner. A short melody will indicate that the PhoneNow feature has automatically switched the hearing instrument to your telephone programme. Initially, you may need to move the telephone receiver slightly to find the best position for reliable PhoneNow activation and good hearing on the telephone.

If you have two hearing instruments with the synchronization function enabled, the volume of hearing instrument on the non-phone ear will be turned down.



Only use recommended cleaning agent to clean the telephone prior to placing the magnet on the phone in order to obtain best possible adherence.



PhoneNow precautions

1. Keep magnets out of reach of pets, children and people with mental disabilities. If a magnet is swallowed, please seek advice from a medical practitioner.
2. The magnet used to supplement the magnetic field of a telephone receiver may affect some medical devices or electronic systems. The manufacturer of any magnetically sensitive devices (e.g. pacemakers) should advise you regarding appropriate safety precautions when using your hearing instrument and magnet in close proximity to the medical device or electronic system in question. If the manufacturer cannot issue a statement, we recommend keeping the magnet or a telephone equipped with the magnet 30 cm (12") away from magnetically sensitive devices (e.g. pacemakers).
3. High distortion during dialing or phoning may mean that the magnet is not in the optimal position relative to the telephone receiver. To avoid the issue, please move the magnet to another place on the telephone receiver.
4. Only use magnets supplied by ReSound.

Telecoil (optional)

A telecoil can be activated by your hearing care professional and accessed through one of the additional programmes. A telecoil picks up a telephone's magnetic signal and converts it to sound. An optional telephone programme may help to improve speech understanding on the telephone. When using a telecoil programme, the receiver of the telephone may need to be held closer to the hearing instrument. The handset of the telephone may need to be moved to slightly different positions in order to find the best reception.

Tele-loop systems

Many places, such as theatres, houses of worship, and schools are equipped with tele-loop systems. When using a telecoil programme with tele-loop systems, sound is picked up directly and may improve speech understanding. If there is no sound from the hearing instruments in a tele-loop system and with a telecoil programme activated, the tele-loop system may not be turned on or is not operating correctly. If a facility is not equipped with a tele-loop system, sitting as close as possible to the front may be helpful.

Direct audio input (optional)

Use of direct audio input (DAI), which enables a direct connection of the hearing instruments to items such as television, radio, and remote microphones, may increase speech understanding for some individuals. The sound source is connected to the hearing instruments by a cable or a wireless FM system to the audio boot. This accessory connects to the bottom of the hearing instruments, and once properly clicked into place, the hearing instruments switch to DAI automatically.






Connecting/Disconnecting audio boots

Connecting audio boots

1. Align the tip of the audio boot with the groove just above the battery compartment and below the model number.
2. Once in place, move the boot in the direction of the battery compartment.
3. Gently click the audio boot onto the hearing instrument.

Disconnecting audio boots

4. Press and hold the button on the front side of the audio boot.
5. Gently remove the audio boot from the hearing instrument.

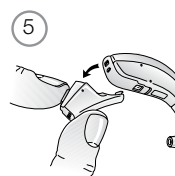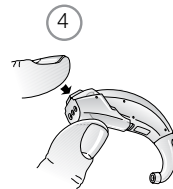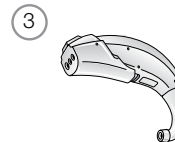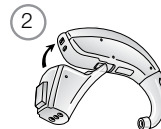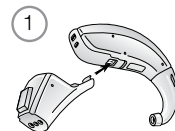



Care and maintenance

Your hearing instrument is protected by a layer of protective, hydrophobic nanocoat material.

Please follow the following instructions to prolong the durability of your hearing instruments:

1. Keep your hearing instrument clean and dry. Wipe the case with a soft cloth or tissue after use to remove grease or moisture. Do not use water or solvents, as these can damage the hearing instrument(s).
2. Never immerse hearing instruments in water or other liquids, as liquids may cause permanent damage to the hearing instruments.
3. Avoid rough handling of hearing instruments or dropping them on hard surfaces or floors.
4. Do not leave hearing instruments in or near direct heat or sunlight, such as in a hot, parked car, as excessive heat can cause damage or deform the casing.
5. Do not wear your instrument while showering, swimming, in heavy rain or in a moist atmosphere such as a steam bath or sauna.
6. If your instrument does get wet, or if it has been exposed to high humidity or perspiration, it should be left to dry out overnight with the battery out and the battery compartment open. It is also a good idea to put the instrument and battery in a sealed container together with a drying agent (desiccator) overnight. Do not use the instrument until it is completely dry. Consult your hearing care professional as to which drying agent to use.
7. Remove your hearing instrument when applying such things as cosmetics, perfume, aftershave, hair spray, and suntan lotion. These might get into the instrument and cause damage.



Daily maintenance

It is important to keep your hearing instrument clean and dry. On a daily basis, clean the hearing instruments using a soft cloth or tissue. In order to avoid damage due to humidity or excessive perspiration, the use of a drying kit is recommended.



Cleaning earmoulds

1. Remove the earmould and attached tubing from the hearing instruments prior to cleaning.
2. Clean the earmould using a mild soap, and rinse with lukewarm water.
3. After cleaning, dry earmoulds thoroughly and remove any residual water and debris from the tubing utilising an air bulb and wire loop.




Note: Earmould tubing may become stiff, brittle, or discoloured over time. Contact your hearing care professional regarding tube changes.



Cleaning thin tubes and domes

1. Remove thin tubes from hearing instruments before cleaning by unscrewing them counter clockwise.
2. Wipe down thin tubes and domes with a damp cloth.
3. In order to clear the thin tube of moisture and debris, push the black cleaning rod through the thin tube, beginning at the end opposite the dome.

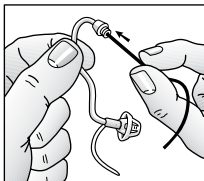


Note: Thin tube and dome systems should be changed every three months or sooner, should the components become stiff or brittle.

How to apply domes

It is recommended that your hearing care professional change domes, as incorrect dome replacement could result in the dome falling out in the ear.

ReSound Standard domes


1. Push the new dome over the flanges on the thin tube.
2. Make sure that the new dome is properly and securely mounted.

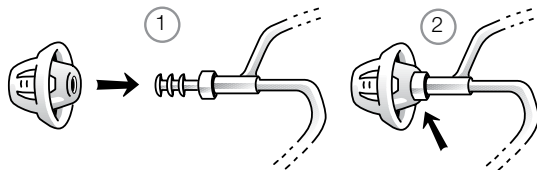

ReSound Tulip domes

Tulip domes are mounted in a similar manner to standard domes, but a few extra steps are required. Tulip domes consist of two “petals”. It is important to note that the largest petal is the outermost petal. To ensure this:

1. Push the largest petal away from the thin tube using a finger. This bends the petal forward.
2. Then push the largest petal back, and it will be placed on top of the smaller petal.

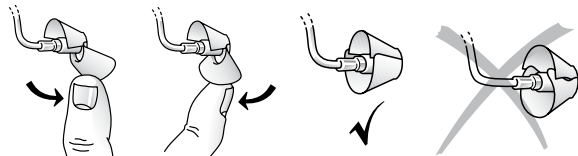

! General precautions

1. Do not leave hearing instruments in the sun, near an open fire, or in a hot, parked car.
2. Do not wear hearing instruments while showering, swimming, in heavy rain, or in a moist atmosphere such as a steam bath or sauna.
3. Should the hearing instrument become moist, remove the battery and place the hearing instrument in a closed container with a drying agent. Your hearing care professional can provide options for drying containers or kits.
4. Remove the hearing instruments when applying items such as cosmetics, perfume, after-shave, hair spray, and suntan lotion.
5. When wireless function is activated, the device uses low-powered digitally coded transmissions in order to communicate with other wireless devices. Although unlikely, nearby electronic devices may be affected. In that case, move the hearing instrument away from the affected electronic device.
6. When using wireless functionality and the devices are affected by electromagnetic interference, move away from the source.
7. Use only original GN ReSound consumables e.g. tubes and domes. Never attempt to modify the shape of the hearing instrument, ear-moulds, or tubing yourself.
8. Only connect ReSound hearing instruments to ReSound accessories intended and qualified to be used with ReSound hearing instruments.

! General warnings

1. Consult a hearing care professional if you think there may be a foreign object in your ear canal, if you experience skin irritation, or if excessive ear wax accumulates with the use of the hearing instrument.
2. Different types of radiation, from e.g. NMR, MRI, or CT scanners, may damage hearing instruments. It is recommended not to wear hearing instruments during these or other similar procedures. Other types of radiation, such as burglar alarms, room surveillance systems, radio equipment, mobile telephones,

contain less energy and will not damage hearing instruments. However, they have the potential to momentarily affect the sound quality or temporarily create strange sounds from hearing instruments.

3. Do not wear hearing instruments in mines, oil fields, or other explosive areas unless those areas are certified for hearing instrument use.
4. Do not allow others to use your hearing instruments. This may cause damage to the hearing instruments or to the hearing of the other individual.
5. Instrument usage by children or mentally challenged persons should be supervised at all times to ensure their safety. The hearing instrument contains small parts that could be swallowed by children. Please be mindful not to leave children unsupervised with this hearing instrument.
6. Hearing instruments should be used only as prescribed by your hearing care professional. Incorrect use may result in hearing loss.
7. Warning to hearing care professionals: Special care should be exercised in selecting and fitting hearing instruments with maximum sound pressure level that exceeds 132dB SPL with an IEC 60711:1981 occluded ear simulator. There may be a risk of impairment of the remaining hearing.
8. Be careful when boarding flights to deactivate the wireless functionality. Turn off your wireless functionality by using the flight mode in areas where radio frequency emission is prohibited.
9. If device is broken, do not use.
10. External devices connected to the electrical input must be safe according to the requirements of IEC 60601-1-1, IEC 60065, or IEC 60950-1, as appropriate (wired connection, f.ex. HI-PRO), SpeedLink).



Note:

* ReSound wireless devices operate in the frequency range of 2.4 GHz - 2.48 GHz.

* ReSound wireless devices include a RF transmitter that operates in the range of 2.4 GHz - 2.48 GHz.

* For use of wireless functionality only use ReSound Unite accessories. For further guidance regarding e.g. pairing, please refer to the user guide of the relevant ReSound Unite accessory.



For ReSound Alera™ TS with Tinnitus Sound Generator

Your ReSound Alera™ TS hearing instrument models also include a Tinnitus Sound Generator function, a tool for generating sounds to be used in tinnitus management programs to relieve suffering from tinnitus. The Tinnitus Sound Generator can generate sounds adjusted to the specific therapeutic needs and your personal preference as determined by your doctor, audiologist, or hearing care professional. Depending on the selected hearing instrument program and the environment you are in, you will sometimes hear the therapeutic sound.

Important Notice to Prospective Users

Good health practice requires that a person with a hearing loss and/or a tinnitus condition have a medical evaluation by a licensed physician (preferably a physician who specializes in diseases of the ear) before using a hearing instrument and/or a sound generator, such as the ReSound Alera™ TS. The purpose of a medical evaluation is to ensure that all medically treatable conditions that may affect hearing and/or tinnitus are identified and treated before the hearing instrument and/or sound generator is used. The sound generator instrument is a tool to generate sounds to be used with appropriate counselling and/or in a tinnitus management programme to relieve patients suffering from tinnitus.

Prescription Use for Sound Generator Users

Please use the device as prescribed by your hearing care professional. Should you develop any side effects from using the instrument, such as dizziness, nausea, headaches, perceived decrease in auditory function or increase in tinnitus perception, you should discontinue use of the device and seek medical evaluation. The target population is primarily the adult population over 18 years of age. This product may also be used with children 5 years of age or older. However, children and physically or mentally challenged users will require training by a hearing care professional or the guardian for the insertion and removal of the device.



CAUTION– Tinnitus Sound Generator

The maximum output of the Tinnitus Sound Generator feature falls into the range that can cause hearing loss according to OSHA regulations.

The user should not use the sound generator for more than eight (8) hours a day when this is set below 90db SPL. Above that level, the device should not be used for more than two hours per day. In no case should the sound generator be worn at uncomfortable levels. Children and physically or mentally challenged users will require guardian supervision while wearing the device.



General warnings/or precautions

Hearing instruments and sound generators can be dangerous if improperly used. Sound generators should be used only as advised by your doctor, audiologist, or hearing care professional. Sound generators are not toys and should be kept out of reach of anyone (especially children and pets) who might cause themselves injury.

Technical specifications **AUDIO SIGNAL TECHNOLOGY**

Digital

AVAILABLE SOUNDS

White noise signal which can be shaped with the following configurations:

- | | |
|-----------------------------|----------------------------|
| - High-pass filter: 500 Hz | - Low-pass filter: 2000 Hz |
| - High-pass filter: 750 Hz | - Low-pass filter: 3000 Hz |
| - High-pass filter: 1000 Hz | - Low-pass filter: 4000 Hz |
| - High-pass filter: 1500 Hz | - Low-pass filter: 5000 Hz |
| - High-pass filter: 2000 Hz | - Low-pass filter: 6000 Hz |

The white noise signal can be modulated in amplitude with an attenuation depth of up to 14dB.



Battery warning information

Batteries, although very small, contain dangerous substances, and should be disposed of carefully. This is for the safety of you and the environment. Please note:

1. Do not attempt to recharge batteries (Zinc Air) which are not specifically designated as rechargeable because they may leak or explode.
2. DO NOT attempt to dispose of batteries by burning them. Used batteries are harmful to the environment. Please dispose of them according to local regulations or return them to your hearing care practitioner.
3. DO NOT place batteries in your mouth. Consult a physician immediately if a battery has been swallowed, as they can be harmful to your health.
4. Keep batteries away from pets, children and individuals who are mentally challenged.
5. Remove the batteries to prevent leakage when the hearing instruments are not in use for an extended period of time.

Troubleshooting Guide

SYMPTOM	CAUSE	POSSIBLE REMEDY
No sound	Not turned on	Turn on by closing the battery door
	Dead battery	Replace battery
	Battery door will not close	Insert battery properly
	Blocked earmould or tube	Clean earmould or tube
	Blocked wax filter	Replace wax filter or consult your hearing care professional
Not loud enough	Incorrect earmould placement	Reinsert earmould
	Blocked earmould or dome	Clean earmould, replace dome, replace filter
	Blocked sound outlet filter	Change filter or consult your hearing care professional
	Change in hearing sensitivity	Consult your hearing care professional
	Excessive ear wax	Consult your physician
	Volume set too low	Increase the volume control if available or consult your hearing care professional

Troubleshooting Guide

SYMPTOM	CAUSE	POSSIBLE REMEDY
Excessive whistling / feed-back	Incorrect earmould placement	Re-insert earmould carefully
	Incorrect dome placement	Re-insert dome
	Excessive ear wax	Consult your hearing care professional
	Feedback control may need adjustment	Consult your hearing care professional
	Earmould tubing worn or damaged	Consult your hearing care professional
	Thin tube connection loose	Change thin tube or consult your hearing care professional
	Hearing instrument settings not optimal	Consult your hearing care professional
Sound distorted / not clear	Weak battery	Replace battery
	Improper earmould or dome fit	Consult your hearing care professional
	Hearing instrument damaged	Consult your hearing care professional
	Hearing instrument settings not optimal	Consult your hearing care professional
Wireless does not work	Possible Root Cause - Device is in flight mode	For devices with push button: Open and close the battery compartment. For devices with-out push button: Open and close the battery door twice within 10 seconds (if Root Cause is device in flight mode)

Technical data

HEARING INSTRUMENT MODEL	MAXIMUM OUTPUT (2cc Coupler / IEC 60118-7)
Models AL977/87-DVIW, AL777/87-DVIW, AL577/87-DVIW, AL477/87-DVIW, AL967-DIW, AL767-DIW, AL567-DIW, AL467-DIW VO988-DW, VO977-DW, VO788-DW, VO777-DW are available in following variants:	
AL967-DIW, AL767-DIW, AL567-DIW, AL467-DIW, ALT967-DIW, ALT767-DIW (Open)	119 dB SPL (typical)
AL977-DVIW, AL777-DVIW, AL577-DVIW, AL477-DVIW, ALT977-DVIW, ALT777-DVIW, ALT577-DVIW (Open)	123 dB SPL (typical)
VO967-DW, VO767-DW (Open)	117 dB SPL (typical)
VO977-DW, VO777-DW (Open)	122 dB SPL (typical)

Technical data

HEARING INSTRUMENT MODEL	MAXIMUM OUTPUT (2cc Coupler / IEC 60118-7)
AL967-DIW, AL767-DIW, AL567-DIW, AL467-DIW, ALT967-DIW, ALT767-DIW	123 dB SPL (typical)
AL977-DVIW, AL777-DVIW, AL577-DVIW, AL477-DVIW, ALT977-DVIW, ALT777-DVIW, ALT577-DVIW	127 dB SPL (typical)
AL987-DVIW, AL787-DVIW, AL587-DVIW, AL487-DVIW, ALT987-DVIW, ALT787-DVIW, ALT587-DVIW	130 dB SPL (typical)
VO967-DW, VO767-DW	122 dB SPL (typical)
VO977-DW, VO777-DW	128 dB SPL (typical)
VO988-DW, VO788-DW	131 dB SPL (typical)

Warranty and repairs

ReSound provides a warranty on hearing instruments in the event of defects in workmanship or materials, as described in applicable warranty documentation. In its service policy, ReSound pledges to secure functionality at least equivalent to the original hearing instrument. As a signatory to the United Nations Global Compact initiative, ReSound is committed to doing this in line with environment-friendly best practices. Hearing instruments therefore, at ReSound's discretion, may be replaced by new products or products manufactured from new or serviceable used parts, or repaired using new or refurbished replacement parts. The warranty period of hearing instruments is designated on your warranty card, which is provided by your hearing care professional.

For hearing instruments that require service, please contact your hearing care professional for assistance. ReSound hearing instruments that malfunction must be repaired by a qualified technician. Do not attempt to open the case of hearing instruments, as this will invalidate the warranty.

Temperature test, transport and storage information

GN ReSound Hearing Instruments are subjected to various tests in temperature and damp heating cycling between -25 C and +70C according to internal and industry standards.

During transport or storage, the temperature should not exceed the limit values of -20C to 60C and relative humidity of 90% RH, non condensing (for limited time). The air pressure between 500 and 1100 hPa is appropriate.

Be aware of information marked with the warning symbol



WARNING points out a situation that could lead to serious injuries,
CAUTION indicates a situation that could lead to minor and moderate injuries.



Advice and tips on how to handle your hearing instrument better.



Equipment includes RF transmitter.



Product is a Type B applied part.



Please ask your local hearing care professional
concerning disposal of your hearing instrument

Worldwide Headquarters

ReSound A/S
 Lautrupbjerg 7
 DK-2750 Ballerup, Denmark
 Tel.: +45 45 75 11 11
 Fax: +45 45 75 11 19
www.resound.com

CVR no. 55082715

United Kingdom

GN ReSound Ltd.
 Kirtlington Business Centre
 Portway
 Kirtlington
 Oxon OX5 3JA
 Tel.: +44 1869 352 800
 Fax: +44 1869 343 466
www.gnresound.co.uk

Australia

GN ReSound Pty. Ltd.
 Unit R1 Regents Park Estate
 391 Park Road
 Regents Park NSW 2143
 Tel.: (free) 1800 658 955
 Fax: +61 2 9743 7472
www.gnresound.com.au

New Zealand

GN ReSound (NZ) Ltd.
 Ground Floor, North Entrance
 4 Fred Thomas Drive
 Takapuna
 Auckland, 0622
 Tel.: (free) 0800 900 126
 Fax: (free) 0800 007 695
www.gnresound.co.nz



Any issues relating to the EU Medical Device Directive 93/42/EEC, or Council Directive 1999/5/EC on Radio Equipment and Telecommunications terminal equipment should be directed to ReSound A/S

ReSound

■■■■■■■■■■
 rediscover hearing